


# Lifestream Field Generator

## Model 100

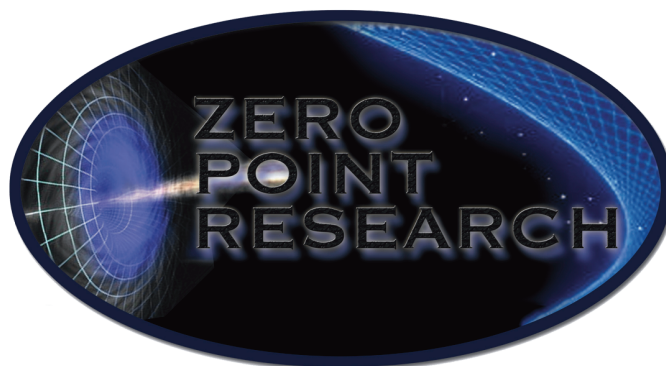

### **!CAUTION!**

This technology is a Radiant Energy Device. It is a new understanding in energy devices. Please take care of this device as it is a highly calibrated instrument and is delicate. Care should be used when using this device. Keep unit clean at all times. Avoid residual build up of elements. Avoid jarring the unit or treating it roughly. Avoid excess weight on surface of the unit. This will cause the unit to fall out of calibration and service will be required.

The *Lifestream Generator* is a combined technological achievement basing research from Nikola Tesla, George Lakhovsky and George W. Van Tassel. There is a combination of over 100 years of research to bring this technology to mankind. The device is based on technology once used during the early 19th Century when electromagnetism was done with EMF (Electro Motive Force) before solid state technology. Faraday called this particular procedure "A conundrum" for the theory of electromagnetism and thus extensive research on our behalf has brought this technology forward to prove its life giving effects to all life on this planet.

The *Lifestream Generator* uses subtle high energy waves to send certain frequencies through the body and in the process raising the cell voltage of all parts of the body thus bringing back a balance to oneself and allowing the body and mind to bring itself back into well being

### **INSTRUCTIONS and PROCEDURES:**

The top of the unit has a built in antenna that allows the energy to radiate out of the unit from one direction.

The ideal position for the experience is 2 to 3 feet from the device using the external antenna and a maximum of 8 people can share the experience at once. You may direct the top of the generator in any direction wanted (as to work on particular areas of the body). During operation the unit will become warm so it is important not to cover the unit with anything and not to run the unit longer than 45 minutes at a time. If continued use is necessary please allow one hour off to 30 minutes on. The generator has an life expectancy of 1200 hours of use. At that time the unit will make unusual noises and start to act erratically. When this occurs please discontinue use and return the unit for service and calibration.

### **DURATION OF USE:**

Depending on age and condition of ones self there is a guide below to recommend the proper duration of the energy experience:

AGE	DURATION	INTERVALS
2 to 15	2 to 5 Minutes	When Necessary
15 to 25	3 to 5	1 per month
25 to 35	5 to 15	1 per week
35 to 45	5 to 15	1 per week
45 to 55	10 to 20	1 per week
55 to 65	15 to 20	1 to 3 per week
65 to 75	15 to 20	1 to 3 per week
75 or Higher	15 to 20	1 to 3 per week

**\* NOTE:** The above reference is for people with normal stress and physical conditions. Additional experiences can be added depending on conditions of the person. Please contact us to obtain a specialized process for your particular needs.